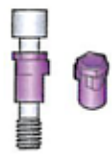






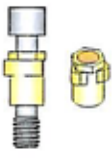






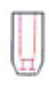






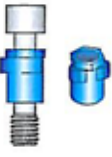

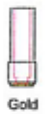










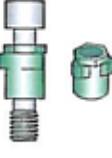


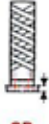










This component is used to convert the internal Tri-nex to an external hex connection. These converters are supplied with a special prosthetic screw and it is essential that this screw is used in conjunction with the prosthesis.

A conventional Tri-nex impression is taken in the chair and a model is cast by the laboratory using the Tri-nex analogue. The converter is then placed in the analogue and the prosthesis can now be manufactured. The prosthesis is supplied on the model to the restorative dentist and the converters are then placed in the same position in the mouth before the prosthetics are introduced.

Converters	Prosthetic Components				Angled Abutments		
<p><b>CTH-35</b> (Hex = 2.54mm)</p>  <p>Packed with *TS-CL-18 Screw</p>	<p>SBN 1 (Hex) SBN 5 (Non-Hex)</p>  <p>Plastic</p> <p>OR</p> <p>GC-EX-34 (Hex) GC-NX-34 (Non-Hex)</p>  <p>Gold</p>	<p>TCBN 1h (Hex) TCBN 1Nh (Non-Hex)</p>  <p>Plastic</p> <p>OR</p> <p>TCBN 5h (Hex) TCBN 5Nh (Non-Hex)</p>  <p>Gold</p>	<p>DN</p>  <p>2 / 3.5 / 5</p>	<p>DBNS12</p>  <p>Scalloped</p>	<p>Flowcharts for Angled Abutment Prosthetics are in the External Hex Connection Product Catalogue CAT-2020</p>		
<p><b>CTH-43</b> (Hex = 2.70mm)</p>  <p>Packed with *TS-CL-20 Screw</p>	<p>SB1 (Hex) SB5 (Non-Hex)</p>  <p>Plastic</p> <p>OR</p> <p>GC-EX-40 (Hex) GC-NX-40 (Non-Hex)</p>  <p>Gold</p>	<p>TCB 1h (Hex) TCB 1Nh (Non-Hex)</p>  <p>Plastic</p> <p>OR</p> <p>TCB 5h (Hex) TCB 5Nh (Non-Hex)</p>  <p>Gold</p>	<p>PKB2H</p>  <p>PEEK</p>	<p>CER-ZR-45</p>  <p>OR</p> <p>CER-ZR-46</p> 	<p>MB</p>  <p>2 / 3.5 / 5</p> <p>OR</p> <p>DB</p>  <p>2 / 3.5 / 5</p>	<p>DBS12</p>  <p>Scalloped</p> <p>OR</p> <p>DBS24</p>  <p>Scalloped</p>	<p>Conical Abutment EB17d OR EB17d-3</p>  <p>GU-CL-9</p> <p>Compact Conical Abutment AMC17d-3 OR AMC30d-4</p>  <p>GU-CL-9</p>
<p><b>CTH-50</b> (Hex = 2.70mm)</p>  <p>Packed with *TS-CL-20 Screw</p>	<p>SBA1 (Hex) SBA5 (Non-Hex)</p>  <p>Plastic</p> <p>OR</p> <p>GC-EX-50 (Hex) GC-NX-50 (Non-Hex)</p>  <p>Gold</p>	<p>TCBA 1h (Hex) TCBA 1Nh (Non-Hex)</p>  <p>Plastic</p> <p>OR</p> <p>TCBA 5h (Hex) TCBA 5Nh (Non-Hex)</p>  <p>Gold</p>	<p>PKBA2H</p>  <p>PEEK</p>	<p>CER-ZR-56</p>  <p>OR</p> <p>CER-ZR-58</p> 	<p>DBA</p>  <p>2 / 3.5 / 5</p>	<p>DBAS12</p>  <p>Scalloped</p> <p>OR</p> <p>DBAS24</p>  <p>Scalloped</p>	<p>Conical Abutment EBA17d OR EBA17d-3</p>  <p>GU-CL-9</p> <p>Compact Conical Abutment ABAMC17D-3 OR ABAMC30D-4</p>  <p>GU-CL-9</p>
<p><b>CTH-60</b> (Hex = 2.70mm)</p>  <p>Packed with *TS-CL-20 Screw</p>	<p>SBBB1 (Hex) SBBB5 (Non-Hex)</p>  <p>Plastic</p> <p>OR</p> <p>GC-EX-60 (Hex) GC-NX-60 (Non-Hex)</p>  <p>Gold</p>	<p>TCBBB1h (Hex) TCBBB1Nh (Non-Hex)</p>  <p>Plastic</p> <p>OR</p> <p>TCBBB5h (Hex) TCBBB5Nh (Non-Hex)</p>  <p>Gold</p>	<p>CER-ZR-67</p>  <p>OR</p> <p>CER-ZR-69</p> 	<p>DBBB</p>  <p>2 / 3.5 / 5</p>	<p>DBBBS12</p>  <p>Scalloped</p> <p>OR</p> <p>DBBBS24</p>  <p>Scalloped</p>	<p>Conical Abutment EBBB17d</p>  <p>GU-CL-9</p> <p>Compact Conical Abutment ABBBMC17D-3 OR ABBBMC30D-4</p>  <p>GU-CL-9</p>	

\* Screws Can Be Purchased Separately

**South Africa**  
Southern Implants (Pty) Ltd.  
Tel: +27 12 667 1046  
Fax: +27 12 667 1029

**Americas / Asia**  
Southern Implants Inc.  
Tel: +1 949 273 8505  
Fax: +1 949 273 8508

**United Kingdom**  
Southern Implants UK  
Tel: +44 208 998 0063  
Fax: +44 208 997 0580

**Australia**  
Henry Schein I Halas Dental  
Tel: +61 2 9697 6288  
Fax: +61 2 9697 6250

**New Zealand**  
Southern Implants Lid NZ  
Tel: 0800 246 752  
Cel: +64 2189 4243

**Germany**  
Southern Implants  
Tel: +49 7121 490 620  
Fax: +49 7121 491 717

**Greece**  
Southern Implants  
Tel: +30 210 898 2881  
Fax: +30 210 595 2543

**Italy**  
S.I.R s.r.l.  
Tel: +39 04 556 8351  
Fax: +39 04 556 8351

**Benelux**  
ProScan bvba  
Tel: +32 11 822 650  
Fax: +32 11 822 651

**Nordic Countries**  
Protera AB  
Tel: +46 31 291078  
Fax: +46 706 150078

**Spain / Portugal**  
Contactodent  
Tel: +351 214 693 332  
Fax: +351 214 693 329

**Turkey**  
Ekodent  
Tel: +212 343 5233  
ellferdal@ekodent.com

**Namibia**  
MF Services  
Tel: +64 61 225 152  
Fax: +64 61 235 630

**Russia**  
Implant Angle Ltd.  
Tel: +7 495 649 1732  
Fax: +7 495 789 9133